

# **Characterization of a Superconducting nanowire single photon detector**

Max Reicherd

Bachelorarbeit in Physik  
angefertigt im Institut für Angewandte Physik

vorgelegt der  
Mathematisch-Naturwissenschaftlichen Fakultät  
der  
Rheinischen Friedrich-Wilhelms-Universität  
Bonn

Juli 2024

I hereby declare that this thesis was formulated by myself and that no sources or tools other than those cited were used.

Bonn, .....  
Date

.....  
Signature

1. Reviewer: Prof. Dr. Sebastian Hofferberth
2. Reviewer: Prof. Dr. Daqing Wang

---

# Acknowledgements

---

I would like to thank ...

You should probably use `\chapter*` for acknowledgements at the beginning of a thesis and `\chapter` for the end.



---

# Contents

---

<b>1</b>	<b>Introduction</b>	<b>1</b>
<b>2</b>	<b>Working principle of SNSPDs</b>	<b>3</b>
<b>3</b>	<b>Faint laser source for detector characterization</b>	<b>7</b>
3.1	Characteristics of faint laser sources . . . . .	7
3.2	Experimental setup . . . . .	8
<b>4</b>	<b>Characterization of our SNSPD</b>	<b>11</b>
4.1	Dark Count Rate . . . . .	11
4.2	Recovery time . . . . .	13
4.3	Efficiency . . . . .	15
4.4	Discussion . . . . .	16
<b>5</b>	<b>Conclusion and Outlook</b>	<b>17</b>
<b>A</b>	<b>Appendix</b>	<b>19</b>
	<b>Bibliography</b>	<b>21</b>



## CHAPTER 1

---

# Introduction

---





## Working principle of SNSPDs

In its essence, the SNSPD consists of four parts, seen in fig. 2.1. A sapphire underground, a detection area made out of a superconducting nanowire in a serpentine winding and placed on the sapphire underground, a gold contact to supply a bias current through the superconducting nanowire and a fiber coupled to the detection area.

The sapphire layer is used to efficiently dissipate the heat when the wire heats up.

Moreover, the set up is first cooled below the critical temperature of the superconductor ( $2 - 3\text{K}$ ), and a bias current is applied to the superconductor that is lower than the critical current.

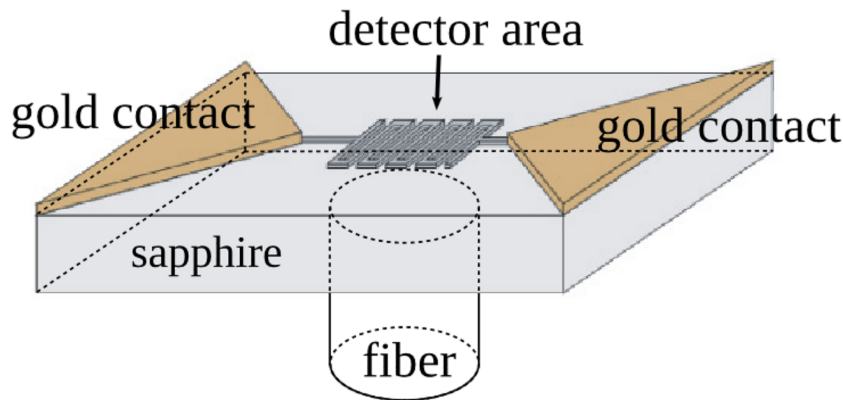


Figure 2.1: (a) Schematic structure of a superconducting nanowire single photon detector [1]

The detection process can be understood along the detection circle in fig. 2.2.

Single photons hitting a superconducting nanowire (ii) and break up individual Cooper pairs. This leads to a local reduction of the critical current below the bias current and in turn to a localised area where the superconductivity is interrupted, this local area forms the so-called "hotspot" (iii). This hotspot forms a resistance area because the bias current exceeds the critical current. In response, the current flows around this hotspot (iv), whereby the local current density in the side areas next to the hotspot again exceeds the critical current, due to a higher current density. This excess also causes a

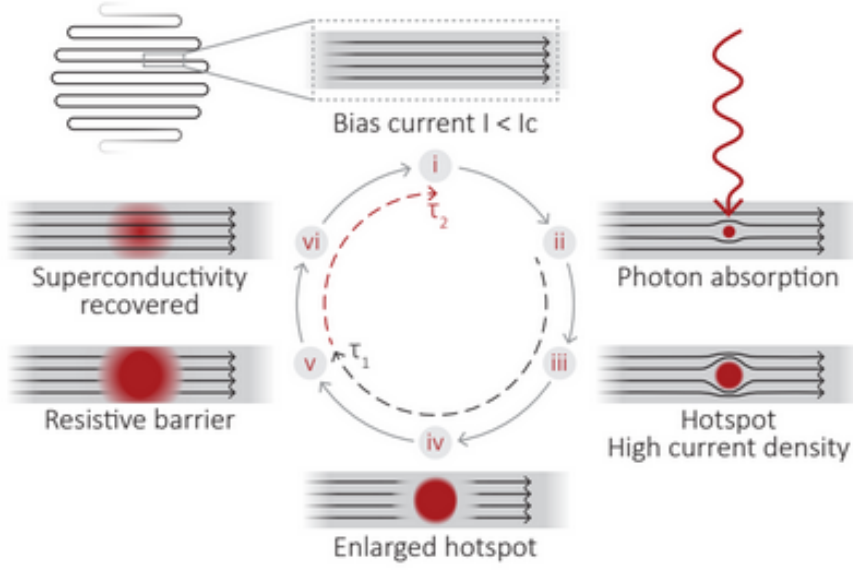


Figure 2.2: Schematic detection cycle of a superconducting nanowire single photon detector [2]

resistance in the side channels due to the critical temperature being exceeded (v). Ultimately, this increase in resistance can be measured in the form of a voltage pulse. The non-superconducting area is then cooled down by the cryogenic environment and returns to the superconducting state (iv—i).

One important technical detail is the fiber coupling of this detector because the efficiency and timing jitter depend on it. Depending on in which polarization the light hits the meander, the efficiency changes. When light hits the wire orthogonally polarized to the wire direction, the photon is less efficiently absorbed, than polarized parallel to the wire.

As seen in fig. 2.3 the coupled fibre in our characterized detector is parallel polarized to the nanowire. Moreover, to cover the whole light, shined out of the fibre, the geometry of the detection zone is constructed as a round plate and has roughly the same diameter as the fibre output (FC/PC). A smaller area would risk not absorbing each photon and a larger area would increase the time leading the signal to the computer unit.

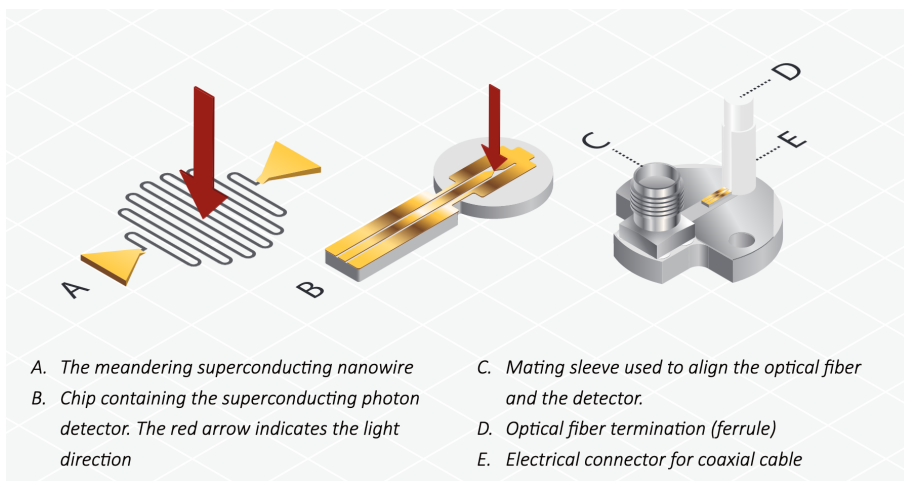


Figure 2.3: Schematic set up of the fiber coupling of a superconducting nanowire single photon detector [Zitieren!](#)



## Faint laser source for detector characterization

Counting photons one has to consider the characteristics of the emitter source as well. Here I focus only on the characteristics of a laser source affecting the characterization I made. For this first, I briefly sum up characteristics of a laser light source and its conditions it gives us for our detector characterization. Further, I make considerations regarding the inherent blindspots of laser light due to the randomness nature of photons in coherent light. Afterward, I introduce the setup I build for characterizing our SNSPD.

### 3.1 Characteristics of faint laser sources

Using a laser source enables us, considering the emitting light as monochromatic beam with angular frequency  $\omega$  and constant Intensity  $I$ . The Photon flux of a laser is defined as the photon average number passing through a cross-section in unit time:

$$\Phi = \frac{IA}{\hbar\omega} = \frac{P}{\hbar\omega} \text{photons } s^{-1} \quad (3.1)$$

where  $I$  is the current of photon,  $A$  the cross-section,  $P$  the laser power and  $\omega$  the angular frequency which depends on the wavelength.

The average number of registered counts  $N(T)$  for a given detection time  $T$  by a detector is given by:

$$N(T) = T\Phi\eta = \frac{PT\eta}{\hbar\omega} \text{photons} \quad (3.2)$$

and hence the registered counts  $\mathcal{R}$  per unit time by:

$$\mathcal{R} = \frac{N(T)}{T} = \eta\Phi = \frac{P\eta}{\hbar\omega} \text{photons } s^{-1} \quad (3.3)$$

where  $\eta$  is the efficiency of the detector system.

This detection count rate is restricted by the largest amount on the dead time of the detector.

$$\mathcal{R}_{\max} \propto \frac{1}{\tau_d} \quad (3.4)$$

Because the used laser has a minimum threshold power which corresponds to a photonrate above the maximum one, one has to attenuate the laser power in order to detect all events.

The photon statistic of coherent light, in our case (in reasonable approximation) of our laser light, is given by poisson statistics. This characteristic stems from discrete nature of photons and hence non-equidistant spacing between photons. Measuring single photons, one has to ensure that we have a neglectable amount of photons in the segment of the deadtime because else our light characteristics inherently forbid us measuring each of the incoming photons.

This is calculated by looking at the propability of measuring one photon per length segment, given by the deadtime. First, we consider one Length segment given by the deadtime  $\tau_d$  and the measurement time  $\tau_m$ :

$$L_d = \frac{c}{\tau_d} \quad (3.5)$$

$$L_m = \frac{c}{\tau_m} \quad (3.6)$$

Through this, we can calculate for a given measurement time the average photon rate per length segment  $L_m$  and further the probability of finding one photon in the particular line segment  $L_d$ :

$$\bar{n} = \Phi \frac{L_m}{c} \quad (3.7)$$

$$p = \frac{\bar{n}}{N} \quad (3.8)$$

Where  $N = \frac{c}{L_d}$  are the subsegments of the measured length segment  $L_m$ .

This enables us to calculate the probability  $p$  of finding  $n$  Photons per deadtime segment  $L_d$  and including this in our measurements.

## 3.2 Experimental setup

Though not measured by ourselves first, it is known from [Zitat](#) deadtime of the detector is around 20-25ns. This gives us a theoretical maximum detection rate of  $\mathcal{R}_{\max} \propto \frac{1}{\tau_d} = 25 - 50\text{MHz}$ . In order to realize the laser attenuation the following setup was build:

The first coupling of the laser light was done in order to operate with the beam on a lower stage, because the laserbeam was due to its construction on an uplifted stage. Afterward the beam passes a pbs to filter the horizontal polarized E field out. Further a galile telescope was build out of one focal and one diffusing lens for minimizing the beam width so it fits fully on the surface of the AOM crystal. The first order of the AOM was set for flexible voltage modulation of the laser. A cover

was used to filter out the first from the zeroth order of the AOM. Then a flip mount was placed, where Neutral density (ND) filters could be placed in and flexibly placed in and out of the laser beam. The ND filters have the function to attenuate the laser light. At the end, before the laser light was again coupled in two waveplates where used to stabilize the light polarization regarding the slow axis of the fibre. Afterwards the light was coupled back into a fibre, so it can be sent to the detector. It was important that the light was coupled in to a APC/PC to FC/PC optical fibre because the detector only had a FC/PC optical fibre input, in order to maintain higher efficiency coupling in the light.

Besides, this optical setup had to be protected from environmental light. For this, the room where the setup was running was shielded with alu foil which has a reflection coefficient of almost 90% at the operating wavelength of 780nm. Moreover a black box was built. It has the function to avoid further environmental light coupling into the fibre. Additionally, the optical fibre running from the optical setup to the detector was shielded with alu foil as well to avoid absorption from the optical fibre.

### Neutral density (ND) filter calibration

We exactly need to know how many photons (related to Power) hit our detector. For this we need to measure the power (or number of photons) of the laser light we shining on the detector. The light we want to shine on the detector is very weak. So weak, we cannot measure it with powermeters or any other measurement devices (except with the uncharacterized detector) This weak power cannot be reached by regulating the laser power down by the laser device (why not? laser minimum power function only with a power higher than required), but has to be done by Neutral Density filters. So our laser power we send to the detector depends on ND filter.

Now there are two challenges relying on ND filters. First, since the fabric values are not precise, the filters need to be calibrated ourselves. Second, one single ND filter cannot attenuate the laser power enough and therefore one has to stack several density filters on each other, since their OD values are adding up.

From these two constraints, we need to calibrate the ND filters carefully, to get precise and accurate OD values of the ND filters.

To get accurate results, we measure the OD value with two methods to rule out systematic errors in our results. The first method measures the OD values, where the ND filters are set in the flip mount as described in graph ???. The second method measures the OD values after the fibre coupling directly in front of the powermeter outside the blackbox.

The attenuation of ND filters is quantized by the Optical density (OD) value. The measurement of the OD values done by sending light on the filter and measuring the power with and without the filter. The OD value is then calculated by the logarithmic value of the inverse transmission value. The transmission value is given by the proportion of the power with the filter compared without the filter.

Afterward, the final OD values of each method were averaged and the corresponding systematic and statistical errors are considered to get the closest to the true OD value for the ND filters.

Based on these OD values each measurement involving the laser was done.

- Basics of photon distribution of laser, attenuation, poisson statistics
- Set up for laser attenuation
- Single Photon Detection paper





---

# Characterization of our SNSPD

---

In literature, four central characteristics have emerged to quantify the quality of single photon detectors and make their performance comparable. These characteristics are the detector efficiency, the dark count rate (DCR), the recovery time and the timing jitter. In this thesis, I focus on the detector efficiency, the dark count rate (DCR) and the recovery time. There are more than these introduced ones, like an after-pulsing, but these are the most important in the context of this thesis.

## 4.1 Dark Count Rate

The dark count rate (DCR) is the rate of measured detection events that were not intentionally sent to the source. It is measured in counts per second and can be caused by statistical fluctuations in the measurement electronics or by scattered or ambient light from the environment. A low DCR is important and allows differentiation between intended detected counts from the source and detected counts by environment and electronic noise. Furthermore, an increase in the count resolution is achieved, enabling the detection of high signal-to-noise ratios at lower input frequencies.

In the context of SNSPDs, the DCR is dependent on the bias current applied to the nanowire. This is due to the fact that if the bias current approaches the critical current, the electronic noise easily exceeds the energy gap required to trigger a pulse. Furthermore, it is important to perform DCR measurements first in the characterisation process. This serves to determine the limit, where general measurements are not distorted by high noise due to high DCR. Consequently, modifications to the DCR configuration serve as the baseline for all subsequent measurements with the detector and the associated setup.

### Measurement setup

In order to evaluate the DCR, it is necessary to perform measurements in two different setups. First, a setup in which the detector is disconnected from the source and covered entirely in darkness. In such a setup, it can be assumed that no photons from the surrounding environment are striking the detector. This allows for the measurement of the DCR triggered by the electronics and the characteristic of the detector's operating mode.

The later means, the probability a signal is triggered by the detector raises if the bias current increases and the difference  $\Delta I = I_c - I_B$  to the critical current decreases. This is due to the lower energy gap its needed exceeding the critical temperature.

The second setup involves connecting the detector to the faint laser source. This is done in order to facilitate a comparison of the results obtained in the two different setups. It is assumed that the first setup has the lowest DCR. A comparison of the two configurations enables the most optimal improvement of the latter towards the former. This ultimately yields the highest signal-to-noise ratio.

## Results and Discussion

Firstly, the DCR is measured with the protection cap in place on the output connector of the detector. This configuration represents the most shielded environment from external light sources and serves as the reference value for the optimal DCR values. The measurement was conducted by sweeping the bias current from 0 to 35  $\mu\text{A}$  in 0.1  $\mu\text{A}$  increments with an integration time of 200 ms.

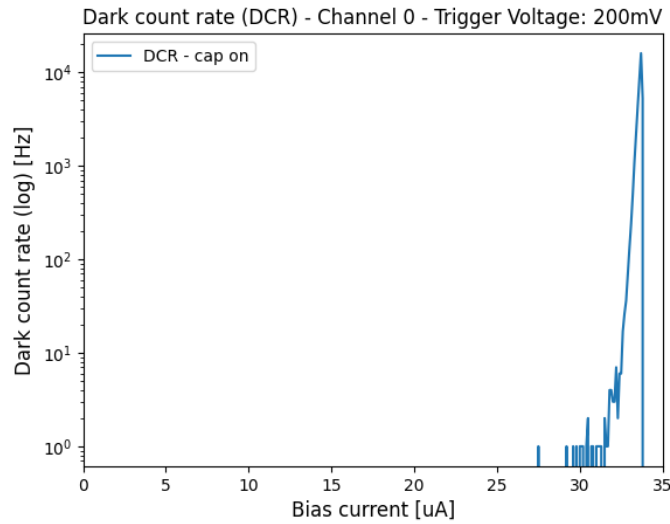


Figure 4.1: Dark count rate (DCR) where the output channel of the detector which was protected by an aluminium cap.

As expected, the DCR raises with decreasing difference  $\Delta I = I_c - I_B$  due to raising probability that weak electronic noises trigger a signal.

Subsequently, the DCR was measured on two further realizations. In the initial configuration, an optical fibre was linked to the detector's output channel and the experiment's output connector. In the second realization, the same connection was maintained, but the setup was fully protected by a self-constructed black box. Furthermore, the optical fibre was enclosed within an aluminium foil layer to prevent photons from coupling through the fibre's cladding.

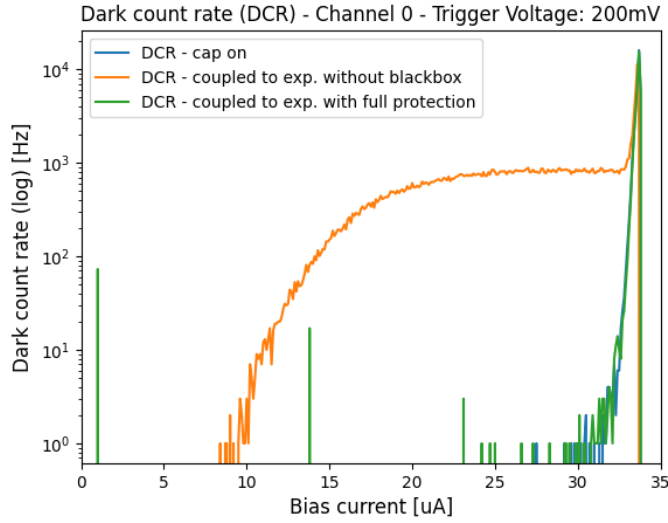


Figure 4.2: Dark count rate (DCR) of both measurement realizations, compared to the DCR with the cap on.

The orange course demonstrates that in the absence of protection, a significant number of photons from the environment are able to couple through various potential pathways, including the fibre cladding and the connection to the experiment.

In contrast, the fully protected setup (shown by the green course) exhibits an analogous course to the setup with the cap on fig 4.1. The DCR of the coupled and protected setup is comparable to the DCR with a cap on. The high peaks at  $\approx 3$  and  $\approx 14\mu\text{A}$  are artifacts resulting probably of still some leakages in the protection. Lastly, all further measurements have to be done at a bias current of  $I_B < \approx 32\mu\text{A}$ . Up to this limit, the DCR is located in the regime of 1-10Hz regime, providing a high lower frequency range for a high signal to a noise ratio.

## 4.2 Recovery time

The concept of the recovery time is visually depicted in fig 4.3. When a photon hits the detector and is absorbed, the efficiency of the detector drops to zero and no further photon can be measured for a certain period of time. This elapsed time is called the dead time  $\tau_d$ . The efficiency then rises again to its original device efficiency  $\eta_G$ . This period is called the reset time  $\tau_r$ . The characteristic curve between the two times forms the start of the increase to full efficiency. The sum  $\tau_{rec} = \tau_r + \tau_d$  of both times forms the recovery time.

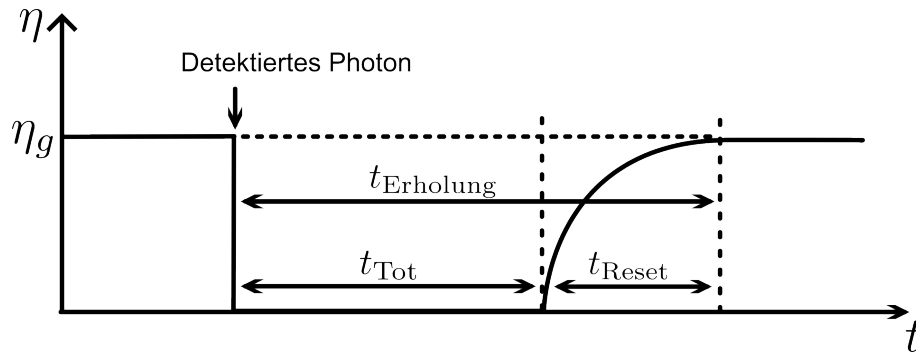


Figure 4.3: Schematic efficiency curve for the detection of a photon[3]. On the Y axis is the efficiency  $\eta$ , where  $\eta_G$  is the device efficiency. On the X axis is the time course of the efficiency

### Measurement set up

The recovery time was measured in the setup seen in fig ?? and measures the recovery time along the same idea provided above. The raw analog signals from the detector were directly transmitted to a time tagger unit (Time Tagger 20) by Swabian instruments (link) with self-adjustable trigger voltages, a device deadtime of 6ns and a maximal counting rate of 9MHz.

This unit enabled the tagging of incoming signals with a time tag, as implied by its name. Subsequently, the tags were used to process the time distances between all signals. The histogram of all distances provides the same efficiency course as previously described, which provides the detectors recovery time for one channel at a determined bias current and trigger voltage.

In order to determine the recovery time of the channel and its optimal configuration regarding bias current and trigger voltage, measurements were conducted for four different bias currents (26, 28,30 and 31.2uA) (below the high noise limit provided by the DCR measurements) and trigger voltages (300, 400, 500 and 600mV).

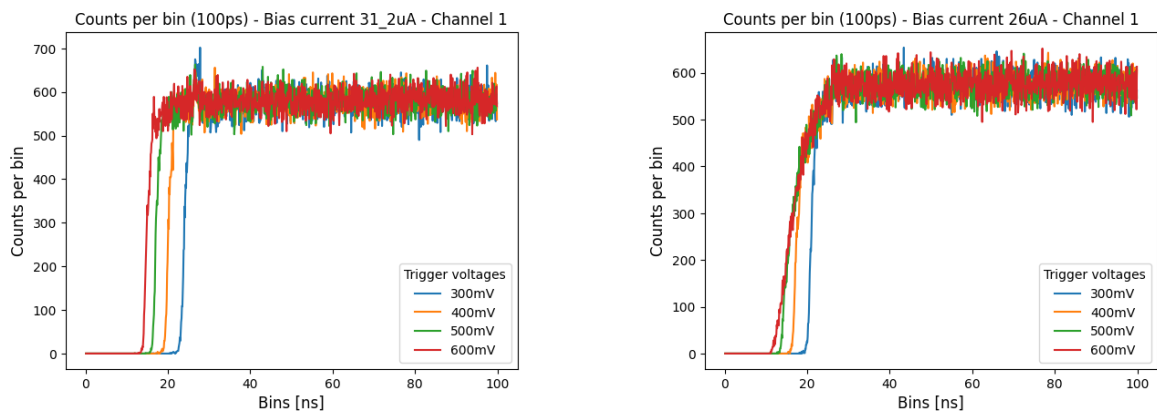


Figure 4.4: isjdias

This was done on the basis that different bias currents influence the recovery of the nanowire back to its superconducting state. Furthermore, the trigger voltage can influence the amount of detected photons. A trigger voltage that is too low will result in the recording of only a single photon, should two photons arrive in quick succession. Like this, the second photon will not be triggered due to the low trigger voltage of the computing device (in this case, the Time Tagger 20). Furthermore, the photon rate was maintained below 9 MHz in order to remain within the operational range of the time tagger.

Additionally, to these measurements analog signals were recorded with an oscilloscope. Reason, for this is to see how the generated pulses look like and whether one can comprehend from them some sequences in our data evaluation. Like this get a better understanding of how the raw output signal changes when parameters as the bias current, the trigger voltage or the input power changes.

The recoverytime was measured with a time tagger connected to the computational unit of the detector. This was done to extract the analog signals and sent them to an alternative to process the data.

### 4.3 Efficiency

In **efficiency**, there are three types of efficiencies that describe independent loss processes in detection. An efficiency can be equated with the probability that a quantum mechanically process under consideration will occur.

These three efficiencies are the coupling efficiency ( $\eta_K$ ), the absorption efficiency ( $\eta_A$ ) and the registration efficiency ( $\eta_R$ ). The graph 4.5 shows schematically where the different loss process in the detection process appear. So, when a photon is sent to a detector via an optical fibre, not all photons can be coupled into the fibre due to material and symmetry properties. The probability of coupling is called the *coupling efficiency*.

When photons hit the detector, there is always a probability that the photon will not be absorbed by the detector. This is described by the *absorption efficiency*.

Finally, there is always a probability that the photon will not be registered by the measuring electronics. This is expressed with the *Registration efficiency*.

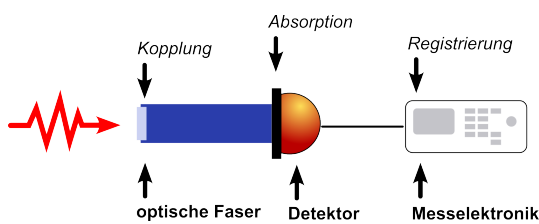


Figure 4.5: Sketch of the components in the detector setup where photon losses appear and consequently a probability ( $\eta_K$ ,  $\eta_A$  or  $\eta_R$ ) has to be considered.

In literature, these terms are brought together in two efficiency terms: the device efficiency ( $\eta_G = \eta_A \cdot \eta_R$ ) and the system efficiency ( $\eta_S = \eta_A \cdot \eta_R \cdot \eta_K$ ). The device's efficiency  $\eta_G$  is that of the device itself and corresponds to photons sent to the detector in a free environment without any fibre

coupling. The system's efficiency  $\eta_S$  also takes into account the coupling losses to the optical fibre. This is the case if the detector is connected to a fibre, as the device properties or the experiment does not allow photon detection in a free environment.

### Measurement set up

The system detection efficiency  $\eta_S$  is measured in different ways, each is pointing out a different variable the efficiency is depending on. Each measurement was done in the setup explained in part 2.2.

First, we wanted to find the polarization axis where the maximum of the light is coupled to the slow axis of the fibre. This is explained by the technical fact that only the slow axis of the fibre is coupled to the output port of the detector. The explanation for this is the maximum  $\eta_A$  explained in part 2.1. By adjusting the laser beam linear with a  $\lambda/4$  plate first and afterwards rotating the  $\lambda$ -half plate in 10 degree steps it was possible to circulate the light axis and hence find the angle configuration where the maximum light was hitting the detector. This is important since measuring the efficiency aligned to a different axis would always put a systematic downshift error on the true efficiency of the detector.

That measurement also confirms Malus law ..

Second, the bias current and trigger voltage dependency was investigated. For this this polarization was aligned to the optimum. Afterwards, the bias current was swept from 0 to 35  $\mu\text{A}$  in 0.1  $\mu\text{A}$  steps and events within 200ms integration time were counted. This was done for four different trigger voltages.

Finally, measurements for different input count rates were done. This was done by exchanging the ND filters combination to get different photon counts. Here I covered up a range from the KHz regime up to the MHz regime. Moreover, this measurement shows the saturation point, where no additional efficiency is acquired by lowering the count rate.

## Results and Discussion

### Results and Discussion

#### 4.4 Discussion

Influence of timing jitter on efficiency and recovery time. - No Afterpulsing - Temperature 2.9 instead of 2.5 K - Photon counters formel für NQO Excelitas SPCS - Note 7

---

## **Conclusion and Outlook**

---





---

## **Appendix**

---

In the appendix you usually include extra information that should be documented in your thesis, but not interrupt the flow.



---

## Bibliography

---

- [1] G. A. Steudle et al.,  
*Measuring the quantum nature of light with a single source and a single detector*,  
Physical review. A, Atomic, molecular, and optical physics/Physical review, A, Atomic,  
molecular, and optical physics **86** (2012),  
URL: <https://doi.org/10.1103/physreva.86.053814> (cit. on p. 3).
- [2] singlequantum.com, *Superconducting Nanowire Single Photon Detectors Operation principle*,  
<https://singlequantum.com/wp-content/uploads/2023/06/Single-Quantum-Operation-Principle-print.pdf>, Accessed: 2023-06, n.d. (Cit. on p. 4).
- [3] K. Shalm et al., *Single-Photon Detector Tutorial*, Computer Scientist QCRYPT, National Institute of Standards and Technology, 2013, URL: <https://qcrypt.github.io/2013.qcrypt.net/contributions/Shalm-slides.pdf> (cit. on p. 14).

